

## bloom U RITE (HP-RECEIVER)

### 4 PERFORMANCE LEVELS 400/300/200/100




The bloom U RITE is based upon the U-Platform with Sound Class Technology for automated classification of sound environments and optimal sound processing determined by the sound classes.

Wireless connectivity via our WidexLink technology, and compatible with the DEX assistive listening devices.

Uses a size 312 battery.

Moderate to severe-to-profound hearing losses.

#### SUGGESTED FITTING RANGE

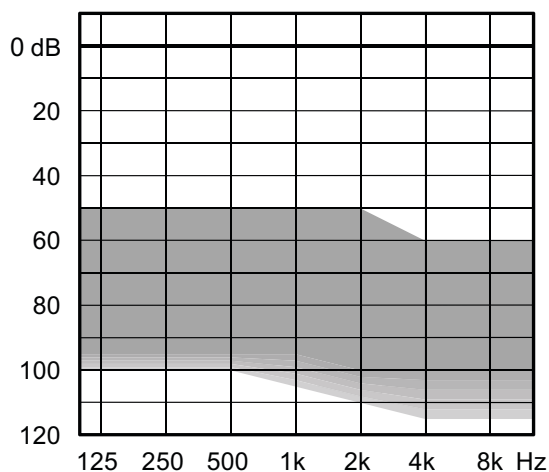

#### STANDARD TECHNOLOGY

- U Platform with Sound Class Technology
- WidexLink wireless connectivity
- InterEar functionality
- Multiple earware options
- Power Saver III technology: Low current consumption

KEY FEATURES	400	300	200	100
Performance	xxxxx	xxxx	xxx	xx
Processing channels	15	10	6	4

#### CONNECTIVITY AND WIRELESS COMPATIBILITY

App for iPhone, Android and new functionality	x Via COM-DEX	x Via COM-DEX	x Via COM-DEX	x Via COM-DEX
WidexLink compatible	x	x	x	x
DEX assistive listening devices (using WidexLink)	CALL-DEX, UNI-DEX, COM-DEX, TV-DEX, RC-DEX, FM+ DEX, PHONE-DEX			

#### UNIQUE FEATURES

FEATURES	400	300	200	100
U platform	•	•	•	•
High-frequency boost	•			
Push button*	U	U	U	U
Wind noise attenuation	U			
Focus Mode	•	•	•	•
Digital Pinna	•	•		
Soft-level noise reduction	•	•	•	•
TruSound Softener	•	•	•	
Speech Enhancer RT	U/IE	IE		
Noise Reduction	•	•	•	•
Sound Class Technology	9 (IE)	5 (IE)	3	1
HD Locator	15	10	6	Broadband
Programs	5	4	3	3
Processing and fine tuning channels	15	10	6	4
ZEN	IE	IE	IE	•
Audibility Extender	U	U	U	U
Variable Speed Compression	U/IE	U/IE	U/IE	U

\* Programmable : Preference Control, program shift or a combination of the two


Protection Class IP58

# BFS U

## MAXIMUM OUTPUT - EAR SIMULATOR

IEC 60118-0


Output dB SPL

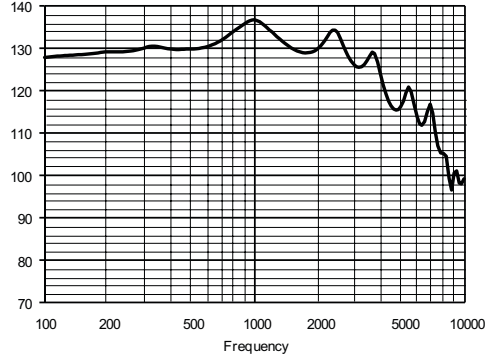

## MAXIMUM OUTPUT - 2CC COUPLER

IEC 60118-7 / ANSI S3.22-2009


Output dB SPL

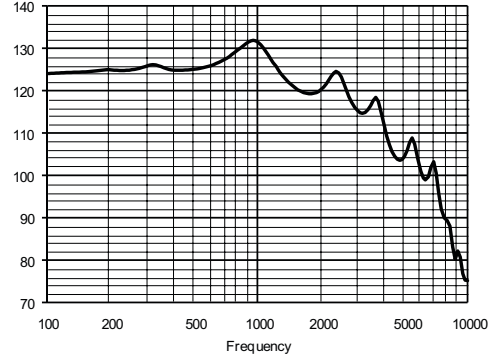

## OUTPUT - EAR SIMULATOR

IEC 60118-0


Output dB SPL

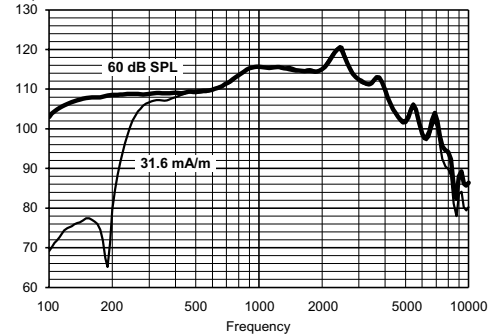

## OUTPUT - 2CC COUPLER

IEC 60118-7 / ANSI S3.22-2009

Output dB SPL

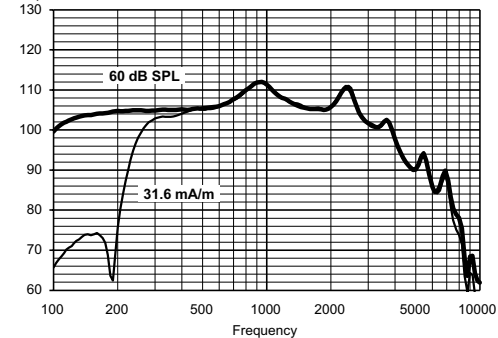

**Technical data** Typical data obtained through standard pure tone measurements. Hearing aid set to Compass Reference Test Gain, unless stated otherwise. Measured using a standard ITE coupler and a custom earmould with Output Extender. For further information, please contact bloom hearing.

		IEC 60118-0	IEC 60118-7 / ANSI S3.22-2009
OSPL90	1600 Hz Peak Average	130 dB SPL 137 dB SPL 132 dB SPL	120 dB SPL 132 dB SPL 126 dB SPL
Acoustic output (Input 60 dB SPL)	1600 Hz Peak Average	115 dB SPL 121 dB SPL 113 dB SPL	105 dB SPL 112 dB SPL 109 dB SPL
Full-on gain (Input 50 dB SPL, Compass Full-on gain)	1600 Hz Peak Average	76 dB 84 dB 75 dB	67 dB 74 dB 70 dB
Telecoil output (Input 31.6 mA/m)	1600 Hz Peak Average	115 dB SPL 121 dB SPL 113 dB SPL	105 dB SPL 112 dB SPL 109 dB SPL
Frequency range		100 Hz - 6000 Hz	100 Hz - 5800 Hz
Harmonic distortion	500 Hz 800 Hz 1600 Hz	<2% <2% <2%	<2% <2% <2%
Equivalent input noise		21 dB SPL	21 dB SPL
Battery drain (stand by)		0.99 mA	0.99 mA
Battery drain		1.06mA	1.10 mA
Battery life / hours (Type 312 Zn-Air, 145 mAh)		135 (>115)	130 (>110)
Mobile phone immunity		IRIL: -28/-7/-9 dB SPL	U-rating: M4/T4