

## bloom U RIC (M-RECEIVER)

### 4 PERFORMANCE LEVELS 400/300/200/100




The bloom U MINI-RIC is based upon the U-Platform with Sound Class Technology for automated classification of sound environments and optimal sound processing determined by the sound classes.

Wireless connectivity via our WidexLink technology, and compatible with the DEX assistive listening devices.

Uses a size 10 battery.

Minimal to severe hearing losses.

#### SUGGESTED FITTING RANGE

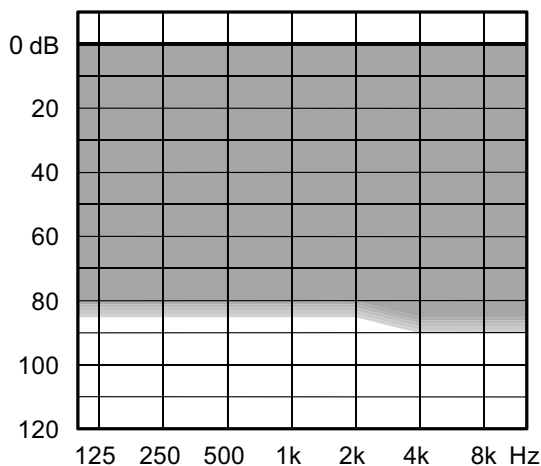

#### STANDARD TECHNOLOGY

- U Platform with Sound Class Technology
- WidexLink wireless connectivity
- InterEar functionality
- Multiple earware options
- Power Saver III technology: Low current consumption

KEY FEATURES	400	300	200	100
Performance	xxxxx	xxxx	xxx	xx
Processing channels	15	10	6	4

#### CONNECTIVITY AND WIRELESS COMPATIBILITY

App for iPhone, Android and new functionality	x Via COM-DEX	x Via COM-DEX	x Via COM-DEX	x Via COM-DEX
WidexLink compatible	x	x	x	x
DEX assistive listening devices (using WidexLink)	CALL-DEX, UNI-DEX, COM-DEX, TV-DEX, RC-DEX, FM+ DEX, PHONE-DEX			

#### UNIQUE FEATURES

FEATURES	400	300	200	100
U platform	•	•	•	•
High-frequency boost	•			
Preference control	U	U	U	U
Wind noise attenuation	U			
Focus Mode	•	•	•	•
Digital Pinna	•	•		
Soft-level noise reduction	•	•	•	•
TruSound Softener	•	•	•	
Speech Enhancer RT	U/IE	IE		
Noise Reduction	•	•	•	•
Sound Class Technology	9 (IE)	5 (IE)	3	1
HD Locator	15	10	6	Broadband
Programs	5	4	3	3
Processing and fine tuning channels	15	10	6	4
ZEN	IE	IE	IE	•
Audibility Extender	U	U	U	U
Variable Speed Compression	U/IE	U/IE	U/IE	U


Protection Class IP58

# BPA U

## MAXIMUM OUTPUT - EAR SIMULATOR

IEC 60118-0


Output dB SPL

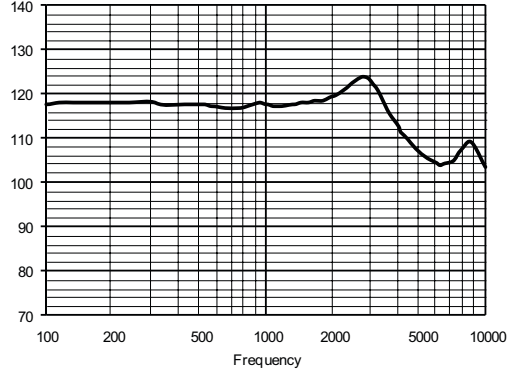

## MAXIMUM OUTPUT - 2CC COUPLER

IEC 60118-7 / ANSI S3.22-2009


Output dB SPL

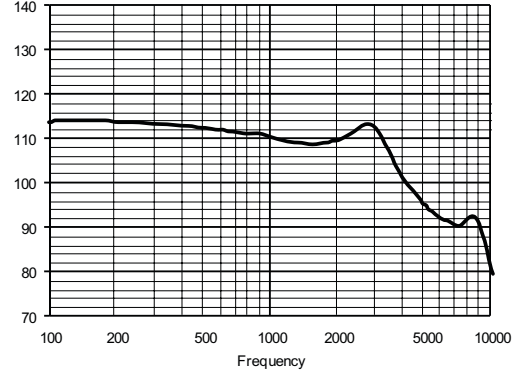

## OUTPUT - EAR SIMULATOR

IEC 60118-0


Output dB SPL

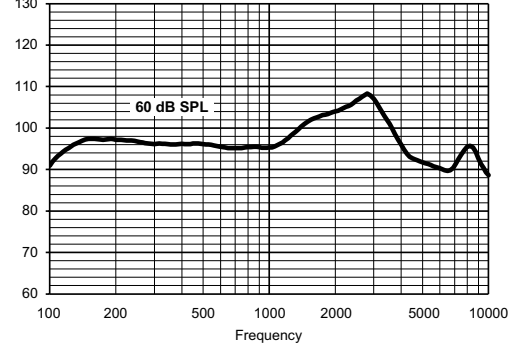

## OUTPUT - 2CC COUPLER

IEC 60118-7 / ANSI S3.22-2009

Output dB SPL

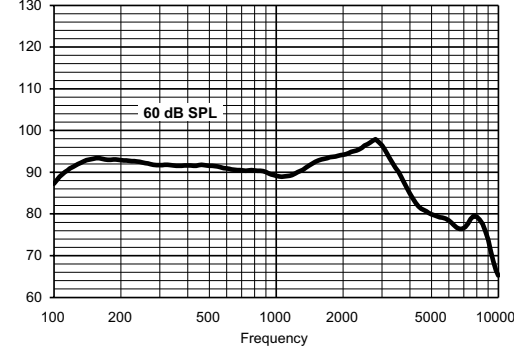

**Technical data** Typical data obtained through standard pure tone measurements. Hearing aid set to Compass Reference Test Gain, unless stated otherwise. Measured using a standard ITE coupler without wax guard. For further information, please contact bloom hearing.

		IEC 60118-0	IEC 60118-7 / ANSI S3.22-2009
OSPL90	1600 Hz	118 dB SPL	109 dB SPL
	Peak	124 dB SPL	114 dB SPL
	Average	118 dB SPL	111 dB SPL
Acoustic output (Input 60 dB SPL)	1600 Hz	102 dB SPL	93 dB SPL
	Peak	108 dB SPL	98 dB SPL
	Average	98 dB SPL	93 dB SPL
Full-on gain (Input 50 dB SPL, Compass Full-on gain)	1600 Hz	63 dB	54 dB
	Peak	69 dB	59 dB
	Average	64 dB	56 dB
Frequency range		100 Hz - 10000 Hz	100 Hz - 9000Hz
Harmonic distortion	500 Hz	<2%	<2%
	800 Hz	<2%	<2%
	1600 Hz	<2%	<2%
Equivalent input noise		22 dB SPL	22 dB SPL
Battery drain (stand by)		0.96 mA	0.96 mA
Battery drain		0.98 mA	1.00mA
Battery life / hours (Type 10 Zn-Air, 90 mAh)		90 (>75)	90 (>75)
Mobile phone immunity		IRIL: -33/-17/-18 dB SPL	U-rating: M4