

GrandAudition RIC with M-receiver

4 Performance levels




The GrandAudition RIC is based upon the E-Platform with a Sound Class controller that handles automatic processing more accurately and faster than before. The GrandAudition RIC use smart technology that learn from the users' preferences and help guide them to a better, more personalised sound. The RIC has an optional ZPower rechargeable solution.

2.4 GHz Bluetooth connectivity with the customisable App allows for direct streaming of audio as well as direct control from smartphones and tablets.

- Direct wireless mobile connectivity (2.4 GHz)
- Multiple wireless connectivity via Apps and DEX assistive listening devices
- Uses an M-receiver
- Uses a size 312 battery
- Protection class IP68 (only non-rechargeable solution)
- Minimal to severe hearing losses

SUGGESTED FITTING RANGE

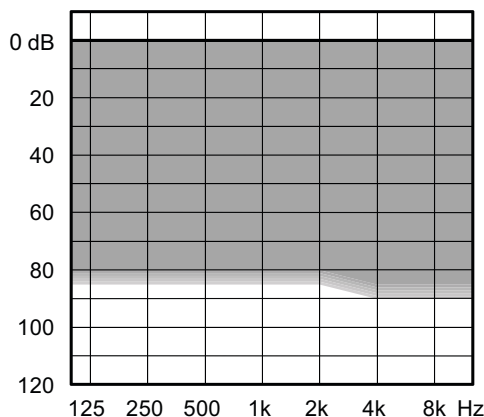

STANDARD TECHNOLOGY

- E-platform with with Sound Class Controller
- Improved open-fit rationales
- Acclimatisation rationales
- Power Saver IV technology: Low current consumption


FEATURES	440	330	220	110
Performance	★★★★★★	★★★★★★	★★★★	★★
Platform	E	E	E	E
SoundSense Adapt	•	•	•	
Adaption manager	•	•	•	•
High-frequency boost	•			
Wind noise reduction	•			
Speech Enhancer RT	RT/IE	IE		
Digital Pinna	•	•		
TruSound Softener	•	•	•	
Soft-level noise reduction	•	•	•	•
Noise Reduction	•	•	•	•
Sound Class Technology 2	11 (IE)	7 (IE)	4	3
HD Locator	•	•	•	
Programs	5	4	3	3
ZEN IE	•	•	•	•
Audibility Extender	•	•	•	•
Preference Control	•	•	•	•
Programmable Push Button*	•	•	•	•
Telecoil	•	•	•	•
ACCESSORIES	440	330	220	110
EVOKE App with SoundSense Learn	•	•	•	
EVOKE App	•	•	•	•
TONELINK App	•	•	•	•
COM-DEX App	•	•	•	•
DEX assistive listening devices**	•	•	•	•
Multiple earware options	•	•	•	•

*Programmable: Preference Control, program shift or a combination of the two
 **Also includes DEX assistive listening devices: CALL-DEX, TV-DEX, COM-DEX, UNI-DEX, RC-DEX, FM+ DEX, PHONE-DEX




*GrandAudition RIC, direct streaming, is compatible with the following devices: iPhone X, iPhone 8, iPhone 8 Plus, iPhone 7, iPhone 7 Plus, iPhone 6, iPhone 6Plus, iPhone 6S, iPhone 6SPlus, iPhone 5S and iPhone SE using iOS 10 or later. Apple, the Apple logo, iPhone, iPad and iPod touch are trademarks of Apple Inc., registered in the US and other countries.


MAXIMUM OUTPUT - EAR SIMULATOR IEC 60118-0

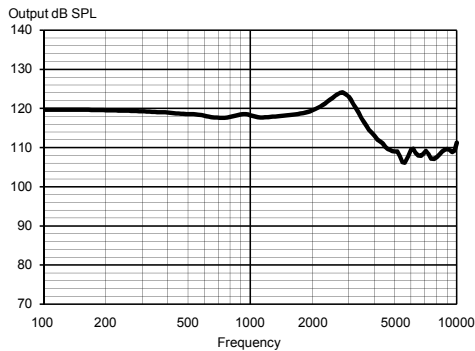


MAXIMUM OUTPUT - 2CC COUPLER IEC 60118-7 / ANSI S3.22-2009

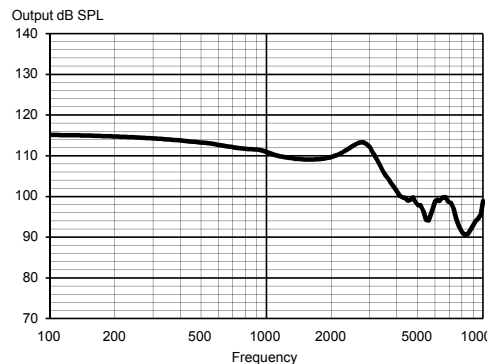

OUTPUT - EAR SIMULATOR IEC 60118-0

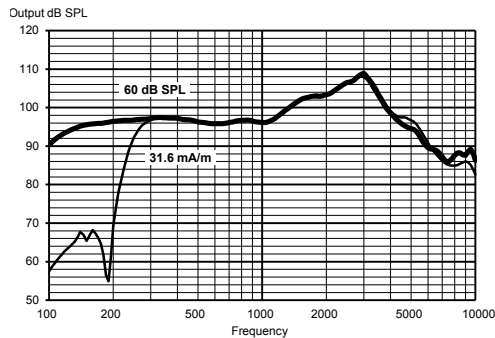

OUTPUT - 2CC COUPLER IEC 60118-7 / ANSI S3.22-2009

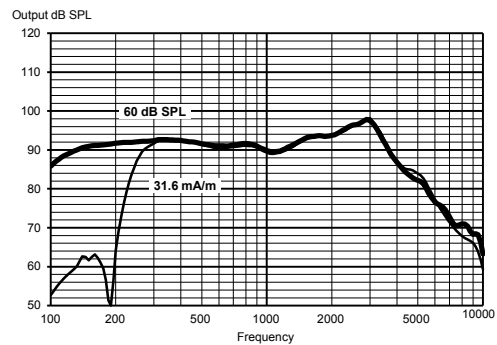

Technical data Typical data obtained through standard pure tone measurements. Hearing aid set to Compass Reference Test Gain, unless stated otherwise. Measured using a standard ITE coupler without wax guard. For further information, please contact GrandAudition.

		EAR SIMULATOR IEC 60118-0:1983 + A1:1994	2CC COUPLER IEC 60118-0:2015 / ANSI S3.22-2014
OSPL90	1600 Hz Peak Average	118 dB SPL 124 dB SPL 118 dB SPL	108 dB SPL 114 dB SPL 110 dB SPL
Acoustic output (Input 60 dB SPL)	1600 Hz Peak Average	103 dB SPL 108 dB SPL 99 dB SPL	93 dB SPL 98 dB SPL 93 dB SPL
Full-on gain (Input 50 dB SPL, Compass Full-on gain)	1600 Hz Peak Average	61 dB 69 dB 62 dB	52 dB 58 dB 55 dB
Telecoil output (Input 31.6 mA/m)	1600 Hz Peak Average	103 dB SPL 108 dB SPL 99 dB SPL	93 dB SPL 98 dB SPL 93 dB SPL
Acoustic frequency range		100 Hz - 10000 Hz	100 Hz - 7000 Hz
Harmonic distortion (typical)	500 Hz 800 Hz 1600 Hz	<2% <2% <2%	<2% <2% <2%
Equivalent input noise		21 dB SPL	21 dB SPL
Battery drain (stand by)		1.00 mA	1.00 mA
Battery drain*		1.02 mA	1.04 mA
Battery life / hours (Type 312 Zn-Air, 170 mAh)* (Type 312 rechargeable, 40 mAh)		165 h 40 h	165 h 40 h
Mobile phone immunity (IEC 60118-13:2016, ANSI C63.19:2011)		IRIL: -21/-19/-9 dB SPL	U-rating: M4/T4

*Battery life in real-life situations depends among other things on the hearing aid features used, streaming time, and the quality of the battery used.